JOCCA


WELLNESS

Ref: 2446

(EN) Instructions manual
(ES) Manual de instrucciones
(FR) Mode d'emploi
(IT) Manuale di istruzioni
(PT) Manual de instruções
(DE) Gebrauchsanweisung
(NL) Gebruiksaanwijzing
(PL) Instrukcja obsługi
(RO) Manual de instrucțiuni
(DK) Brugsanvisning

DESCRIPTION / DESCRIPCIÓN / DESCRIPTION / DESCRIZIONE /
DESCRIÇÃO / BESCHREIBUNG / BESCHRIJVING / OPIS PRODUKTU /
DESCRIERE / BESKRIVELSE

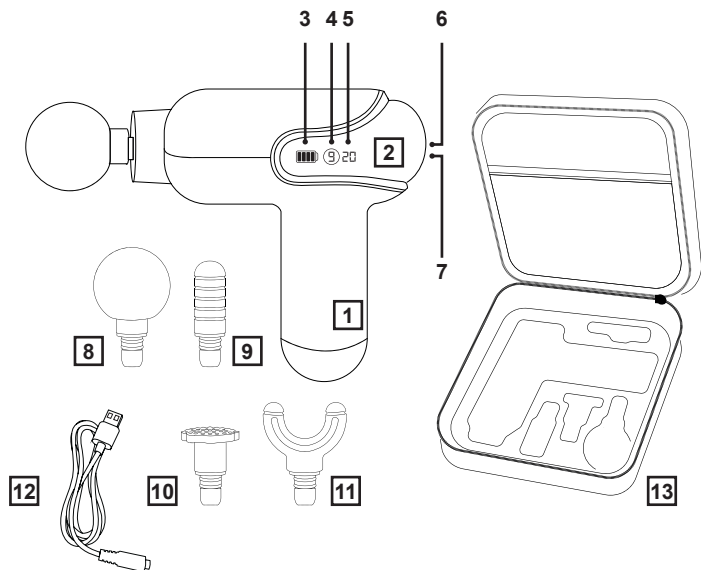EN/DESCRIPTION

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MASSAGE GUN

SAFETY WARNINGS

- These instructions are highly important, please read this instruction manual before using the device and save it for future reference.
- JOCCA is not responsible for any consequences resulting from improper use of the product or any other use not described in this manual.
- Before using the product for the first time, unpack it and check that it is in perfect condition. If not, do not use the product and contact JOCCA, as the warranty covers any damage from origin or manufacturing defects.
- This product is intended for domestic use only. Do not use the device outdoors. Do not expose it to weather conditions such as rain, snow, sun, etc. that may impair its operation.
- This device can be used by children aged from 8 years and above, and people with reduced physical, sensory or mental capabilities or lack of experience or knowledge if they are supervised or instructed in the appropriate and safe use of the appliance and understand the hazards involved. Children should not play with the appliance. Cleaning and product maintenance shall not be carried out by children, unless they are at least 8 years old and under supervision.
- Keep the product and the power cord away from children under 8 years old.
- Children should be supervised to ensure that they do not play with the product.
- Should you need to put the product down while using it, always do so with the device switched off and on a stable, horizontal and smooth surface.
- Before connecting the device to the mains supply,

check that the cord and plug are in perfect condition. Also, check that the voltage indicated on the product is compatible with your electrical installation.

- If your product has been dropped or water or other liquid has been spilled on it, do not use or attempt to repair it. Contact JOCCA's technical service.

- If you detect that the power cord, plug or adaptor is damaged or soaked with water or other liquid, do not use the appliance, or attempt to replace the cord, plug or adaptor. Contact directly authorised JOCCA technical service, which will replace it correctly in order to avoid any danger.

- Do not immerse the product or any part of it in water or any other liquid.

- Always use JOCCA accessories and spare parts for your safety and the best possible performance of the appliance.

- Do not use the appliance with wet hands or feet, or with bare feet.

- Do not leave the appliance running unattended.

- Consult your doctor before using the device if you are in any of the following situations:

- If you are pregnant or have just given birth

- If you suffer from heart disease or have a pacemaker.

- If you suffer from fractures, osteoporosis, or any muscle, spine or joint ailment.

- If you suffer from a malignant tumour.

- Use the massage gun on clean, dry skin. It is recommended to use it over your clothes, pressing gently and moving it over the skin, about 60 seconds for each area.

- Use the massage gun only on the soft tissues of your body without pain or discomfort. Avoid parts above the

neck, with less muscle coverage and bones.

- Do not use the device if you have bruises or if you feel pain or discomfort.
- Do not obstruct or place objects in the ventilation holes.
- Do not bring fingers, hair or other body parts close to the shaft or the back of the massage head to avoid being pinched.

BATTERY WARNINGS

- Remember that batteries must be deposited in the specialized containers for battery collection.
- If you want to throw out the device, always remove the batteries before doing it.
- Never place the battery compartment near fire or a heat source.
- Do not dispose of batteries in fire or place the control unit in the sun.
- If the battery is leaking, never handle it directly with your hands.
- The battery must be removed from the appliance for safe disposal.

INSTRUCTIONS FOR USE

- Before using the device for the first time, fully charge the battery.
- Place the desired head into the massage gun.
- Press and hold the ON/OFF button for 3 seconds to turn the appliance on.
- The appliance has 9 speeds. In addition, if you exert a little pressure on the area to be massaged, the gun will work more intensely.
- To change the speed, press the ON/OFF button until the desired speed is selected. The countdown indicator will show the remaining usage time.

the speed at which the massage gun is operating.


- To switch off the appliance, press and hold the ON/OFF button for 3 seconds.
- The appliance is protected against overheating, so it will automatically switch off after 20 minutes of continuous use. The countdown indicator will show the remaining usage time.

Massage heads and functions:

The massage gun includes 4 heads with different shapes and functions:

	BALL HEAD For large muscle groups such as the back, legs, buttocks, pectorals...
	BULLET HEAD To access deep muscle points or relieve tension in areas such as the palm of the hand or plantar area.
	FLAT HEAD To relax and tone large muscle areas.
	FORK HEAD For skirting and massaging bony areas such as the neck area, around the spine and the Achilles tendon.

- Before changing the heads, make sure that the massage gun is switched off.
- To change the heads, simply pull the head away from the main body and insert the desired head into the slot as shown below.

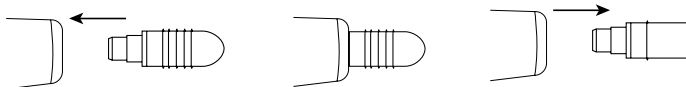Charging

- To charge the massage gun battery, connect one end of the cable to the charging port and the other to a 5V DC USB adapter (not included). The battery indicator will light up on the screen.
- The battery is fully charged when the four lines of the battery indicator are lit.
- Do not use the massage gun while it is charging. Always disconnect it from the mains before use.

CLEANING AND STORAGE

- Before cleaning, make sure that the appliance is turned off and disconnected from the mains supply.
- Clean the outside of the appliance with a damp cloth or sponge and then with a dry cloth.
- Do not use abrasive detergents or products that could damage the appliance.
- NEVER IMMERSE THE APPLIANCE AND THE POWER CORD IN WATER OR ANY OTHER LIQUID.
- Store the device in a dry and safe place.

TECHNICAL DATA

Charge: 5V $\equiv$ 1A (by USB)

7.4V $\equiv$ (Li-Ion)

Lithium battery 1500 mAh

This product complies with the directives 2014/30/EU EMC, 2011/65/EU and 2015/863/EU ROHS.

PRODUCT RECYCLING – ENVIROMENT

Do not dispose of this product at the end of its useful life. Take it to a designated waste collection point (subject to the 2012/19/EU WEEE directive).

HOW TO REMOVE THE BATTERIES:

If the product contains any batteries, they must be removed before disposing of the product by following the steps below:

- Ensure that the batteries are completely discharged before removing them from the appliance.
- The appliance must be switched off and disconnected from the mains when the battery is removed.
- Open the casing of the appliance, previously removing the fixing screws with the help of a screwdriver.
- Identify the battery and cut its wires (red and black). To avoid the risk of short-circuit with the possible residual charge of the battery, first cut the red wire and protect the conductive part with insulating tape. Then cut the black wire in the same way.
- Safely remove the battery from its location.
- Remember that batteries must be disposed of in special, authorised containers.

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