MASS FOOT

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INFORMATIONS TECHNIQUES

Tension / fréquence: 220-240V~, 50/60Hz

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EN Please read its instruction carefully before using the massage, especially safety attentions, so that you can use it correctly.

Besides, please keep the instruction after you read it in case you need refer to it someday.

Thanks for your cooperation!

Safety Attentions

• **This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.**

- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children unless they are aged from 8 years and above and supervised.
- Keep the appliance and the cable out of reach of children.
- Do not try to repair the appliance on your own. Always contact an authorized technician. To avoid the exposure to danger, always have faulty cable be replaced only by the manufacturer, by its customer service or by qualified person and with a cable of the same type.
- The devices are not intended to be operated by an external timer or a remote control system.

1. Persons who are not allowed to use the massage

 WARNINGS	• Pregnant women, patients with heart attack. • Patients whose foot has skin disease.
--	--

2. Something you have to know when you are using the massage.

 ATTENTIONS	• Children are not allowed to use the massage unless with the adults' help. • Persons who are not sensitive to heat can't use the massage for too long every time. • Don't stuff anything to the open parts. • Don't pull the massage on the ground in case it damage the massage. • Don't put your center of gravity on one feet. • The massage is intended for persons' life and it can't be used in other ways.
--	---

 FORBIDDEN	• Stop using the massage if you find it damaged. • Don't stand on the massage or put something heavy on it. • Stop using the massage if it is abnormal.
 NO DEMOUNTING	• Don't take apart the massage by yourself, or it may break down. • Please contact with our after-sale service department if you have any questions.
 PULL OUT	• Don't pull out the attaching plug with wet hands. • Pull out the attaching plug after you using it. • Pull out the attaching plug when there is a power cut.

3. Concerning the attaching plug and power line


 ATTENTIONS	• Don't use the power supply whose voltage is higher than rated voltage, or it may cause fire. • Pull out the attaching plug instead of power lines. • Do use the dedicated power lines to replace the damaged one.
--	---

4. Operating environment


 ATTENTIONS	• Don't put the massage in a wet or dusty place. • Don't put the massage near somewhere in a high temperature. • Don't use the massage when it thunders.
--	--

Functions

1. Sketch map

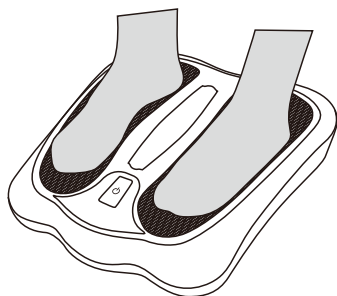2. The name of parts

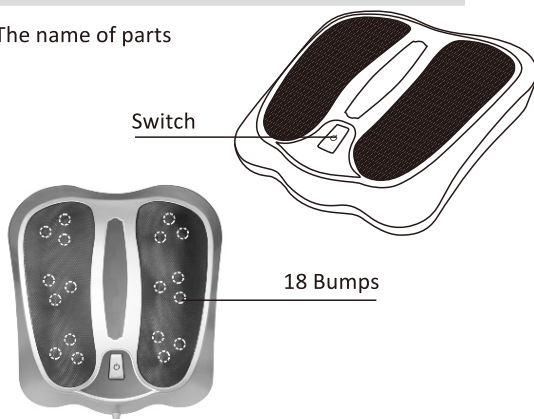3. Features

- Deep kneading and shiatsu massage
It can stimulate foot reflection zone to relieve the pain and relax the channels and collaterals.
- New massage head with 18 rotatable bumps.
You can enjoy the massage as if a professional massagist does it.
- Thermal therapy with far infrared It helps to promote cell metabolism.

Usages

Operation

1. Insert the attaching plug before using. (alternating voltage : 220V Frequency:50Hz)
2. Sit down and put your foot on the massage.
3. Press the switch once and you can enjoy the kneading massage.(The blue light will shine.)
4. Press the switch twice and you can enjoy the kneading massage with heating. (The red light will shine.)
5. Pressing the switch three times means turning off the massage.

! **Attentions**

- Don't use the massage over 20 minutes. (It don't have timing devices.)
- You'd better use the massage with socks on foot.
- Heating cannot be used alone.

Maintenance

1. Cleaning

- Before cleaning, pull out the attaching plug and wait until it become cold.
- Please clean the crust of the massage with a dry rag.
- Don't clean the massage in any liquid.
- Don't use brushes, petrol, cleanser or diluent during cleaning.

2. Storage

- Put the massage at a safe and airy place if you won't use it for a long time.
- Don't hang the massage with wires.

Trouble Shooting

If you meet problems when you are using the massage, please try to resolve it according to the tips below.

Abnormal situations	Possible reasons	Solutions
It doesn't work.	1. It can't energize. 2. You don't turn on the switch.	1. Insert the attaching plug again. 2. Turn on the switch.
Stop working suddenly.	The electricity is cut off.	Check whether the power interface is in order.

If you still cannot resolve the problem, please pull out the attaching plug and contact with our after-sale service department.

Technical Parameters

Voltage Rating: 220-240V~

Frequency Rating: 50/60Hz

Power Rating: 31W

Keep it in case you need it.

It may cause some pollution if you don't deal with the waste correctly.

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