
designed by **BRAUN**

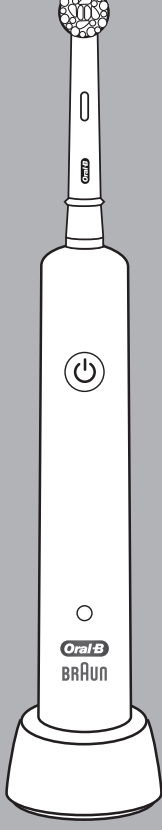Charger / Chargeur	Type 3757
Handles / Manches	Type 3710
	Type 3708
	Type 3791

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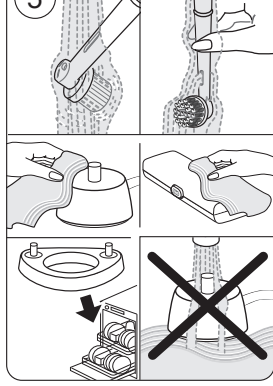
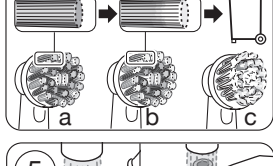
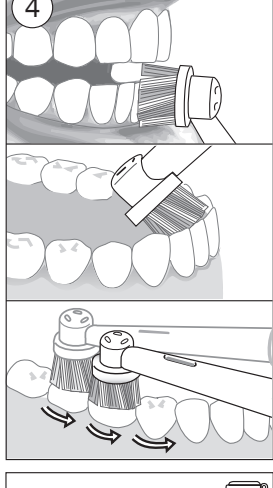
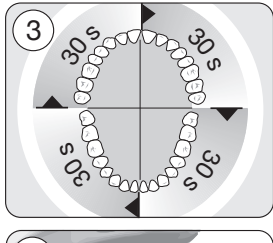
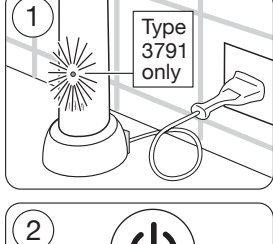
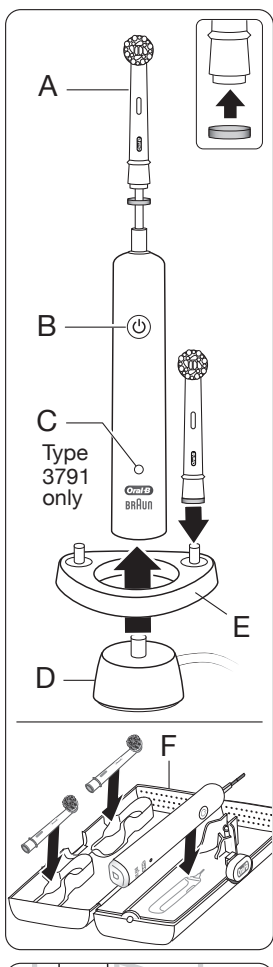
---

## Welcome to Oral-B!

Before operating this toothbrush, please read these instructions and save this manual for future reference.

# IMPORTANT

1. Periodically check the entire product/cord/accessories for damage. A damaged or non-functioning unit should no longer be used. If the product/accessories is damaged, take it to an Oral-B Service Center.  
The supply cord cannot be replaced. If the cord is damaged the appliance should be discarded. Do not modify or repair the product. This may cause fire, electric shock or injury.
2. Never insert any object into any opening of the appliance.
3. Usage by children under age 3 is not recommended. Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
4. For tooth brushing only. Do not chew or bite on brush. Do not use as a toy or without adult supervision.
5. Cleaning and maintenance shall not be performed by children.
6. Children shall not play with the appliance.
7. Use this product only for its intended use as described in this manual. Do not use attachments (e.g. brush head refills) which are not recommended by the manufacturer.
8. Only use charger provided with the appliance.

# WARNING

1. Do not place the charger/charging stand, power supply or charging travel case in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
2. This appliance contains batteries that are non-replaceable. The appliance has no user serviceable parts. Do not open and disassemble device except when disposing of the battery. Follow handling and disposal instructions of the battery/device described in this manual.
3. When unplugging, always pull the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
4. If stickers are provided, follow instructions on the sticker sheet when applying on handle. Ensure valve on back of the handle is not covered.
5. This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution. For hygienic reasons, use individual brush handles and brush heads for each family member.
6. Small parts may come off, keep out of reach of children.
7. To avoid brush head breakage which may generate a small parts choking hazard or damage teeth:
  - Prior to each use, ensure the brush head fits well. Discontinue use of the toothbrush, if brush head no longer fits properly. Never use without a brush head.
  - If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible.
  - Replace brush head every 3 months or sooner if brush head becomes worn.
  - Clean brush head properly after each use (see «Cleaning recommendations» paragraph). Proper cleaning ensures the safe use and functional life of the toothbrush.

---

## Description

- A Brush head
- B On/off switch (mode selection)
- C charge/low-charge indicator light only for Type 3791
- D Charging unit

## Accessories (depending on model):

- E Brush head holder
- F Travel/Travel Charging case (design may vary)

**Note:** Content may vary based on the model purchased.

**To identify your purchased type please refer to the bottom/back of your toothbrush.**

## Specifications

For voltage specifications please refer to the bottom of the charging unit.

## Charging and Operating

This toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

Plug the charging unit (D) into an electrical outlet and place the toothbrush handle on the charging unit.

Full charge depends on Model:

### Type 3710:

A full charge takes about 16 hours and enables up to 8 days of regular brushing (twice a day, 2 minutes)

### Type 3708 and 3791:

A full charge takes about 22 hours and enables up to 10 days of regular brushing (twice a day, 2 minutes)

### Type 3791 only:

The green charge indicator light (c) flashes while the handle is being charged. Once it is fully charged, the light turns off.

Note: After deep discharge, the green light might not flash immediately; it can take up to 10–15 minutes.

If the rechargeable battery is low the red charge indicator light (c) is flashing for some seconds when turning your toothbrush on/off.

## Operating

For everyday use, the handle can be stored on the plugged-in charging unit to maintain it at full power. Overcharging is impossible. To maintain the maximum capacity of the rechargeable battery unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

## Using your toothbrush/Brushing technique

Wet brush head and apply any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching on the appliance (picture 2).

When brushing your teeth with one of the **Oral-B oscillating-rotating brush heads** guide the brush head slowly from tooth to tooth, spending a few seconds on each tooth surface.

With any brush head start brushing the outsides, then the insides and finally the chewing surfaces. Brush all four quadrants of your mouth equally. Do not press too hard, simply let the brush do all the work (picture 4). You may also consult your dentist or dental hygienist about the right technique for you.

In the first few days of using any electric toothbrush, gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist.

## Brushing modes (depending on model)

### Type 3710:

«**Sensitive**» – Gentle, yet thorough cleaning for sensitive areas

«**Daily Clean**» – Standard for daily mouth cleaning

### Type 3708:

«**SensitivePlus**» – Super gentle cleaning for extra sensitive areas

«**Sensitive**» – Gentle, yet thorough cleaning for sensitive areas

«**Daily Clean**» – Standard for daily mouth cleaning

### Type 3791 only:

«**Daily Clean**» – Standard mode for daily cleaning

«**Sensitive**» – Gentle, yet thorough cleaning for sensitive areas

«**Whitening**» – Polishing for occasional or everyday use

«**Gum Care**» – Gentle massage of gums

## Pressure Sensor

### Type 3791 only

To promote optimal brushing your toothbrush has a pressure control feature installed. If too much pressure is applied, the movement of the brush head will continue but its pulsation will stop. In addition, you will also hear a different sound while brushing.

Note: The sequence and availability of brushing modes may vary depending on model purchased.

To switch to another mode, press the on/off button (B). Switch off the toothbrush by holding the on/off button until the motor stops.

## Timer (depending on model)

An incorporated timer memorizes the elapsed brushing time (depending on model) (picture 3).


«**2-Minute**» timer: A short stuttering sound indicates the end of the professionally recommended 2-minute brushing time.

«**Professional**» timer: A short stuttering sound at 30-second intervals reminds you to brush all four quadrants of your mouth equally. A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time.

## Brush heads

Oral-B offers you a variety of different Oral-B brush heads that fit your Oral-B toothbrush handle.

### Oral-B CrossAction brush head

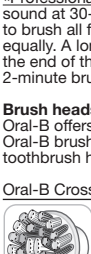Features perfectly angled bristles to remove plaque, even in hard-to-reach areas.

### Oral-B Gum Care brush head

Features a combination of regular and ultra-thin bristles for advanced cleaning and gum protection.

### Oral-B Sensitive Gum Care brush head

Features extra soft bristles for gentle cleaning of sensitive areas.

Most Oral-B brush heads feature INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the color will fade approximately within 3 months, indicating the need to exchange your brush head. If the bristles splay before the color recedes, you may be exerting too much pressure on teeth and gums (picture 4).

We do not recommend using the Oral-B «FlossAction» or Oral-B «3D White» brush head with braces. You may use the Oral-B «Ortho» brush head, specifically designed to clean around braces and wires.


## Cleaning recommendations

After use, rinse the brush head for several seconds under running water with the handle switched on. Switch off and remove the brush head. Rinse both parts separately under running water and wipe them dry before reassembling the toothbrush. Unplug the charger before cleaning. Charging unit (D) and travel case (F) should be cleaned with a damp cloth only. Never place charging unit in water. Store the clean and dry toothbrush/brush heads in the travel case. Brush head holder (E) is dishwasher-safe (picture 5).

Subject to change without notice.

## Environmental notice

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your local region.

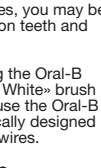
## Battery removal at the end of the product's useful life


Do not disassemble the product except when disposing of the battery. Run down battery before removing. When taking out the battery for disposal, use protective

clothing (gloves, safety glasses, etc.) to avoid direct skin contact with batteries that may be subject to leakage. Use caution not to short the positive (+) and negative (–) terminals. Insulate terminals (e.g. with tape) and put the battery in a plastic bag for disposal.

### Type 3710 and 3708
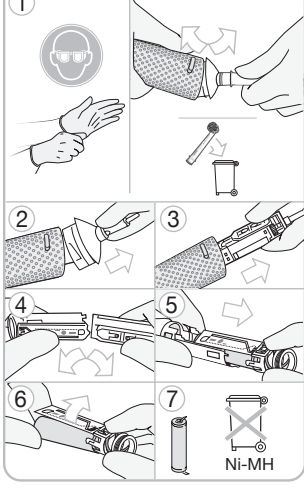
### Type 3791

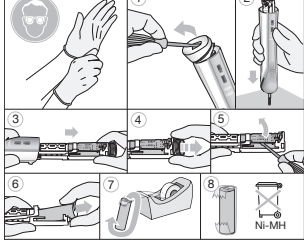Open the handle as shown, remove the battery and dispose of it according to local environmental regulations.

Caution! Opening the handle will destroy the appliance and invalidate the warranty.

### Satisfaction Guaranteed or Your Money Back:

Try it with no risk for 60 days!

If you are not 100% satisfied, call 1-877-593-4199 to receive a full refund. Return the product in its original packaging along with the original cash receipt postmarked within 60 days of purchase to receive a prepaid card in the amount of purchase. Terms of prepaid card apply.

Your Oral-B toothbrush must be returned in its original packaging and be postmarked no later than 60 days after the date of purchase. Please allow 6-8 weeks for handling. In the U.S., the shipment cost to return the product is the responsibility of the consumer. Sorry: No refunds will be made for any Oral-B product damaged by accident, neglect, unreasonable use or lost in the post.

Terms and conditions apply. For further details please contact our customer service.

### Warranty replacement brush heads

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the use of non-Oral-B replacement brush heads.

Oral-B does not recommend the use of non-Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high-quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

## FOR USA ONLY

### Limited 2-Year Warranty

Oral-B warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail.

If the appliance exhibits such a defect, Oral-B will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Oral-B Service Center. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Oral-B. This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL ORAL-B BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential

damages, so the above limitation may not apply to you. This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

### How to obtain service in the USA

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

*To obtain service:*

- A. Carry the product in to the authorized Oral-B Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

*To obtain original Oral-B replacement parts and accessories:*

Contact an authorized Oral-B Service Center.

Please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

*To order original Oral-B replacement parts and accessories (including replacement brush heads) by phone:*

Please call 1-800-566-7252.

*To speak to an Oral-B Consumer Service Representative:*

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